


SomnoDent® is a mandibular advancement device (MAD) created on the basis of the patient's bite registration and consists of upper and lower dental plates.

The **b - flex Comfort Liner** maximizes patient comfort. The resulting high levels of patient satisfaction and acceptance contribute substantially to the therapeutic success of the treatment.

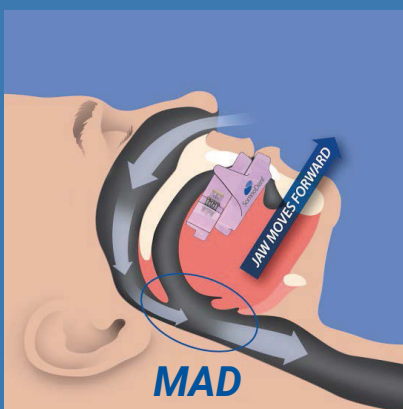

The patented mechanism of our MAD gently moves the lower jaw into the desired therapeutic position and keeps the airway open during sleeping

SomnoDent®

MAD vs CPAP

| | | |
|------------------------------|-------------------------------------|-------------------------------------|
| Permits normal mouth opening | <input checked="" type="checkbox"/> | <input type="checkbox"/> |
| Allows speech and drinking | <input checked="" type="checkbox"/> | <input type="checkbox"/> |
| Easy to travel with | <input checked="" type="checkbox"/> | <input type="checkbox"/> |
| Noiseless | <input checked="" type="checkbox"/> | <input type="checkbox"/> |
| Inconvenient face mask | <input type="checkbox"/> | <input checked="" type="checkbox"/> |
| Limited sleep positions | <input type="checkbox"/> | <input checked="" type="checkbox"/> |
| Complicated cleaning | <input type="checkbox"/> | <input checked="" type="checkbox"/> |
| Requires power connection | <input type="checkbox"/> | <input checked="" type="checkbox"/> |



www.somnomed.com

SomnoMed products are registered in:
EU, UK, Switzerland, United States, Canada,
Australia, South Korea, Japan, Hong Kong

SomnoMed products are patented
and patient matched medical devices.

SomnoMed AG is ISO-13485 certified.

SomnoMed AG

CH-6300 Zug
Switzerland




SOMNOMED'S INNOVATIVE PRODUCT RANGE


www.somnomed.com


Somnomed offers patients, physicians and sleep dentists an evidence-based, comfortable alternative to CPAP for the first time. Patients worldwide have been treated with a SomnoDent device – returning their dreams and the energy to fulfill them. SomnoDent products are indicated for the treatment of mild to moderate obstructive sleep apnea (OSA) and for patients that are intolerant for other therapy options.

Our devices are clinically proven in treating mild to moderate obstructive sleep apnea, supported by medical evidence, innovative in design and provide a comfortable fit for patients.

Intraoral Scanning

 Scan and cloud-send patient's digital impression and protrusive bite directly

 No messy or costly impression materials

 More comfortable for patients - no gagging

 Faster product turn-around time

 Ensures precision and predictability in production

SomnoDent Features: More Comfortable: SomnoDent devices are custom fitted and more comfortable for patients and makes the therapy efficient due to its high tolerance/compliance.

High Quality: We use the highest quality acrylic which won't discolor, and offer a 2- year warranty on our devices. *Please see details at our website.

Clinically Proven: SomnoMed stands by our quality promise by manufacturing devices that are proven and validated by more than 17 clinical studies.

Device Selection: We offer a wide range of sleep apnea devices for different patient types.

Patient Awareness: We inform patients about OSA solutions.

Compliance: Customers rate SomnoDent products as a 4 on a scale of 1-5.

Practice Support: Our sales team is actively connecting sleep dentists with sleep physicians, increasing referrals.

Technical Support: Experienced Lab and Customer Service staff for all technical and product choice inquires.

Marketing Support: Patient Materials, website & advertising support.

SomnoDent® CLASSIC



Classic offers a hard acrylic surface and easy retention adjustment. Provides a high level of durability and retention. Simple to adjust chairside.

SomnoDent® FLEX



For the ultimate in comfort and flexibility Flex offers the standard in a comfortable oral device with a one of a kind soft b-flex Comfort Liner that is heat cured that facilitates bonding between layers and a reliable titration system.

Provides a comfortable "first time fit" for patients with minimal adjustment needed.

Superior device retention even for patients with short teeth, little undercut or crown, bridge work.

SomnoDent® FUSION



For the ultimate innovative comfort and flexibility, a longer titration range and fast titration Fusion offers premium comfort and convenience.

Combination of a removable wing system for quick titration and the basic fine calibration screw technology.

Usable instantly with b-flex Comfort Liner material for a consistent retention and quality.

SomnoDent® EDENTULOUS



Solution for partial and complete edentulous.

The b-flex covers the full maxillary palate and is retained by posterior palatal seal.

Occludes with a lower splint using a two piece coupling mechanism.

Range of motion is preserved without jeopardizing the retention of the upper splint with a downward pull.

A minimum of 4-5 teeth or 2 implants in the lower arch are necessary.

SomnoDent® AVANT



SomnoDent Avant is digitally designed for dentists who want a slim, durable device with a quality fit and finish. Our CAD-CAM manufacturing system ensures a precise and consistent fit. With an innovative titration system and our proprietary b-flex Comfort Liner, the AVANT is our slimmest most comfortable device, particularly suitable for a patient with a narrow arch.

Durable SomnoDent Avant Advancement Straps provide simple calibration in 1mm increments

An initial set of 10 straps provided with your case with additional sizes available in small, medium, large (10x per set) to cover every patient's complete advancement range

Completely digital process produces a precisely fit device with exceptional strength

Backed by SomnoMed's dependable 2-year warranty

Morning Repositioner



Hands-free bite trainer.

Provides a patient-specific, comfortable and compressible raised splint that directs clenching forces straight down.

Uniquely designed for patients to perform isometric exercises that should minimize the possibility of mandibular positioning forward.

SomGauge

For quickly and easily capturing the most accurate bite registration, SomGauge is the first gauge of its kind to allow for an increased measurement of the vertical along with an adjustable ramp.



Provides amplified impression accuracy and device integrity and comfort for the patient.

SomTabs



For gentle and effective cleaning of SomnoDent devices without harming the appliance or producing a bad after taste.

Exclusively manufactured for SomnoDent devices. Does not contain Sodium Perborate